

# Compat ella®

## Compat Ella® Universal Set

The Compat Ella® Universal Set is intended for administration of enteral nutrition, fluids or medication and is suitable for use with the Compat Ella® feeding pump.

### Features at a glance

The Compat Ella® Universal Set features an **ENPlus** nutrition connector, a **universal screw cap** with air filter, a **3-way** medication port and an **ENFit®** distal end (patient side tube connection). A separate **bottle hanger** is also included.


**Built-in Safety:** All Compat Ella® administration sets feature patented fail-safe **Auto-Anti-Free-Flow (AAFF)** technology. **AAFF** is incorporated into the cassette and is always automatically activated to ensure patient safety.

**Safer Connectors:** All Compat Ella® administration sets include **ENFit®** connectors made according to **ISO 80369-3**. The **ENFit®** and **ENPlus** connectors are designed to be **IV Incompatible** to promote **patient safety** by helping to avoid accidental misconnections.

**Food Safe Materials:** All Compat Ella® administration sets are **DEHP, BPA and Latex-Free\***.

**Compatibility:**

- Nutrition containers: ENPlus, glass and plastic wide-neck containers, crown cork glass bottles
- Enteral feeding tubes: ENFit®
- Medication dispensers: ENFit®

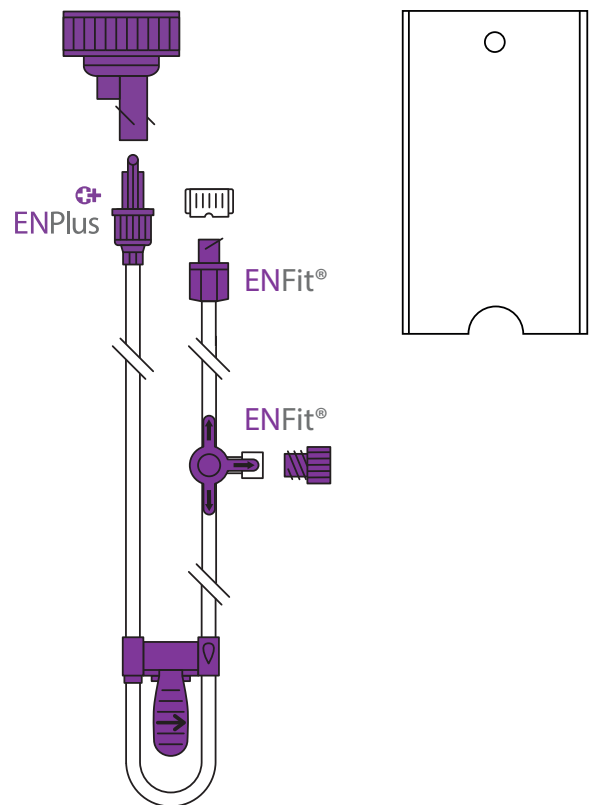


### Ordering Information

Product Name	Product Code	Qty/Box
Compat Ella® Universal Set	12298206	30

Have you met Compat Ella®?

# Compat Ella® Universal Set

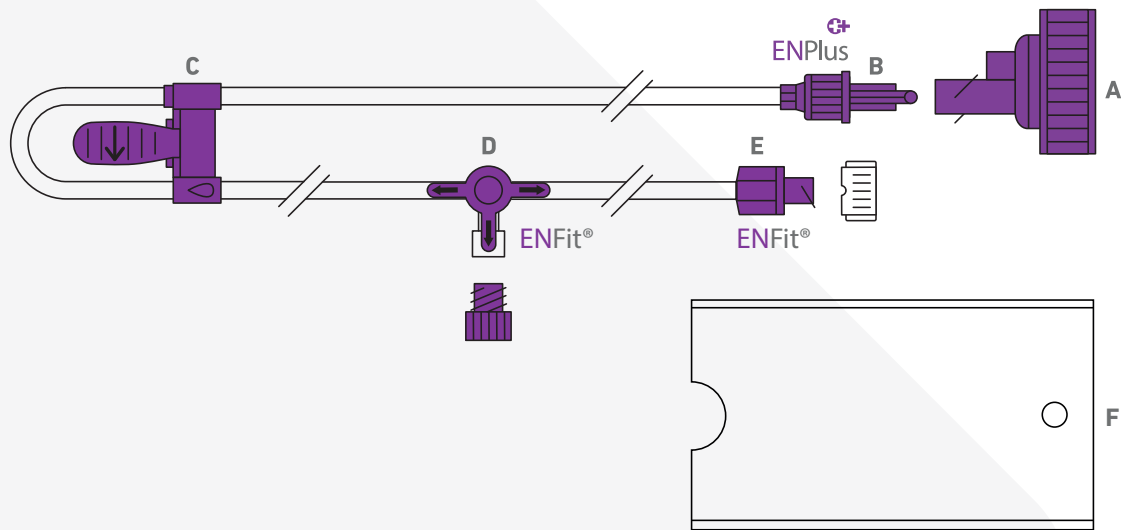

## Components

- A** Universal screw cap with air filter and ENPlus port for use with glass and plastic wide-neck containers (DIN 168-1) or crown cork glass bottles (EN 14634 and EN 14635)
- B** ENPlus 'enteral only' proximal end (nutrition connection) with revolving lock
- C** Compat Ella® pump cassette with integrated and patented Auto-Anti-Free-Flow (AAFF) safety feature
- D** ENFit® 3-way medication port with cap
- E** ENFit® distal end with cap (patient side tube connection)
- F** 500 mL bottle holder LDPE (Low Density Polyethylene)

## Specifications

<b>Sterile:</b>	Sterilized with Ethylene Oxide.
<b>Packing:</b>	Individually packed in cases of 30 pieces (39.5 x 29.5 x 14 cm).
<b>Storage:</b>	Store at room temperature. Keep away from sunlight and humidity.
<b>Shelf life:</b>	3 years from production date. Expiry date on the package.
<b>Duration of use:</b>	24 hours. Single Use.
<b>Tube Length:</b>	Distance from nutrition connection to patient side tube connection: 230 cm. Distance from nutrition connector to pump: 50 cm.
<b>Tube diameters:</b>	Nominal inner diameter: 3.0 mm. Nominal outer diameter: 4.15 mm.
<b>CE Mark:</b>	Complies with the European Medical Device Directive 93/42 EEC. A CE Certification has been issued by an accredited Notified Body: TÜV SÜD (CE0123).
<b>Ecological Footprint:</b>	The 'Grüne Punkt' symbol signals that Nestlé Health Science assumes its producer responsibility and contributes to the financing of the recycling.
<b>Materials:</b>	Tubing: PVC and silicone. All materials DEHP, BPA and Latex-free*. Additional details on request.

To obtain further information please contact your Nestlé Health Science representative or go to [www.compatella.com](http://www.compatella.com)

\*The parts do not exceed the maximum concentrations of 0.1% total DEHP (Di (2-Ethyl Hexyl) Phthalate) content. None was intentionally added during the manufacturing process.

\* Unless stated otherwise, all trademarks are owned by Société des Produits Nestlé S.A or used with permission. Legal Manufacturer: Cedic Srl. Made in Italy.

NHSM00406 REV.1 12DEC2017



# Compat ella®